

PART OF THE **DURAN**® BOTTLE SYSTEM

CLICK AND CARRY

DURAN® GL 45 BOTTLE
CARRYING SYSTEM



**DURAN
WHEATON
KIMBLE**

Excellence in your hands

DURAN® GL 45 BOTTLE CARRYING SYSTEM

The new modular Bottle Carrying System fits standard laboratory bottles that have a GL 45 neck finish and is compatible with four bottle sizes: 2, 5, 10 and 20 litres. It is ergonomically designed, and offers health & safety benefits to the user and allows easy and safe transport of bottles around the laboratory or production facility. It makes pouring liquids from the larger bottles more comfortable and safer. All materials are fully autoclavable.

Carrying Handle

- Carrying handle is suitable for bare or gloved hands
- Handle strap is flexible enough to allow pouring

Tightening Strap

- With indicator arrow
- End modified to prevent unthreading

Pouring Handle

- Ergonomically shaped
- Anti-slip grip



Silicone Support Base

- Structural support for base of larger 5, 10 and 20 litre bottles
- Protection of glass bottle base from impact during lifting or handling
- Acts as a spacer to protect glass from contact with Metal Dolly
- Moulded-in holes minimize liquid traps and allow maximum steam penetration during autoclaving



GL 45 Bottle Collar

- The collar is securely attached by pushing down over the GL 45 thread with 4 teeth that lock under transfer ring.
- The collar is only removable by cutting – so attachment is 100% secure in use
- A cut line indicator for removal of collar when required

Horizontal Straps

- Durable stitched construction
- Made from "anti-slip" material




Lot Specific Retrace Code

for traceability and certification



Optional Extra: Metal Dolly

- Stainless steel dolly for 10 or 20 litre bottles (2 sizes only)
- Fully compatible with carrying system
- Easy to clean, crevice-free construction
- Optional handle is available for easier handling

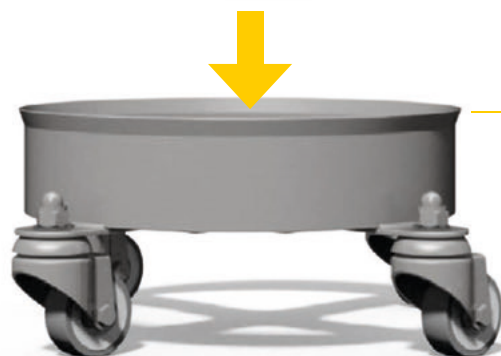

FOR EASY AND SAFE CARRYING & POURING



Safe and controlled pouring using handle

- Safer pouring of liquids using a second lower "Pouring Handle" located on the side strap (Available on the 5 litre carriers only)



Easy and safe carrying using handle

- Handle attaches at the neck and fully supports the base

Fast and secure tightening

of the polypropylene straps with two side strap adjusters for the larger 5, 10 and 20 litre bottles



2 L bottle carrier in 4 colors

for color coding



4 sizes available:
2 L / 5 L / 10 L / 20 L

Find out more:
www.duran-bottle-system.com





DWK Life Sciences GmbH
Hattenbergstraße 10
55122 Mainz
Germany

Phone: +49 6131 - 1445 4131
Fax: +49 6131 - 1445 4016
sales@DWK-LifeSciences.com
www.DWK-LifeSciences.com

