

Minivette® POCT

Collection and dispensing of small blood samples
for Point-of-Care tests




- Easy sample recovery due to optimal capillary features
- Precise dispensing of small whole blood volumes
- Prevents spillage during transfer to test card
- Available in 10 µl, 20 µl, 30 µl, 50 µl, 100 µl and 200 µl volumes
- Preparations: Neutral, Lithium-Heparin and K₃ EDTA




Pre-analytics

Minivette® POCT


Description	Color Code	Packaging	Order No.
Minivette® POCT, 10 µl, neutral	□	200 pcs./inner box / 2.000 pcs./case	17.2111.010
Minivette® POCT, 20 µl, neutral	□	200 pcs./inner box / 2.000 pcs./case	17.2111.020
Minivette® POCT, 30 µl, neutral	□	200 pcs./inner box / 2.000 pcs./case	17.2111.130 (upon request)
Minivette® POCT, 50 µl, neutral	□	200 pcs./inner box / 2.000 pcs./case	17.2111.050
Minivette® POCT, 100 µl, neutral	□	200 pcs./inner box / 2.000 pcs./case	17.2111.100
Minivette® POCT, 200 µl, neutral	□	150 pcs./inner box / 1.500 pcs./case	17.2111.200
Minivette® POCT, 20 µl, Lithium-Heparin	■	200 pcs./inner box / 2.000 pcs./case	17.2112.020
Minivette® POCT, 20 µl, Lithium-Heparin	■	200 pcs./inner box / 2.000 pcs./case	17.2112.120
Minivette® POCT, 30 µl, Lithium-Heparin	■	200 pcs./inner box / 2.000 pcs./case	17.2112.130 (upon request)
Minivette® POCT, 50 µl, Lithium-Heparin	■	200 pcs./inner box / 2.000 pcs./case	17.2112.050
Minivette® POCT, 50 µl, Lithium-Heparin	■	200 pcs./inner box / 2.000 pcs./case	17.2112.150
Minivette® POCT, 100 µl, Lithium-Heparin	■	200 pcs./inner box / 2.000 pcs./case	17.2112.100
Minivette® POCT, 100 µl, Lithium-Heparin	■	200 pcs./inner box / 2.000 pcs./case	17.2112.101
Minivette® POCT, 200 µl, Lithium-Heparin	■	150 pcs./inner box / 1.500 pcs./case	17.2112.200
Minivette® POCT, 20 µl, K ₃ EDTA	■	200 pcs./inner box / 2.000 pcs./case	17.2113.020
Minivette® POCT, 20 µl, K ₃ EDTA	■	200 pcs./inner box / 2.000 pcs./case	17.2113.120
Minivette® POCT, 30 µl, K ₃ EDTA	■	200 pcs./inner box / 2.000 pcs./case	17.2113.130 (upon request)
Minivette® POCT, 50 µl, K ₃ EDTA	■	200 pcs./inner box / 2.000 pcs./case	17.2113.050
Minivette® POCT, 50 µl, K ₃ EDTA	■	200 pcs./inner box / 2.000 pcs./case	17.2113.150
Minivette® POCT, 100 µl, K ₃ EDTA	■	200 pcs./inner box / 2.000 pcs./case	17.2113.101
Minivette® POCT, 200 µl, K ₃ EDTA	■	150 pcs./inner box / 1.500 pcs./case	17.2113.201

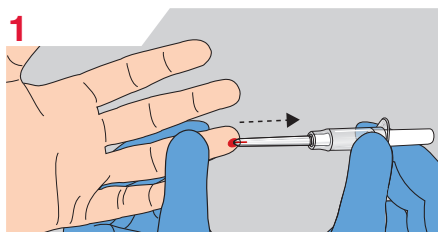

Hold the Minivette® POCT in a horizontal, slightly declined position, under the wings. The air ventilation hole at the end of the piston should not be covered, nor should the piston be pushed down. Collect the blood drops using the capillary tip ensuring there are no air bubbles present.

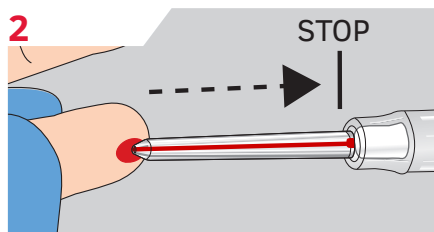


Blood collection is complete when the blood has reached the filter at the other end of the capillary tip.

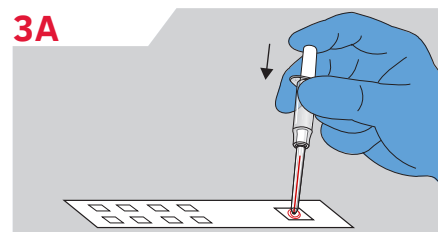

Place the capillary tip on the test field and slightly press down the piston to dispense the whole blood sample.

All patient blood samples must be treated as potentially infectious.

Wear gloves!

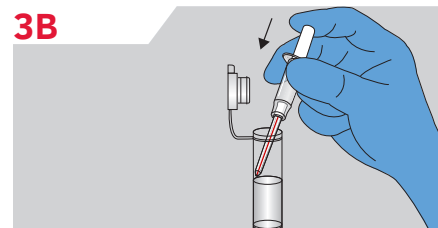

Alternatively, dispense the sample into a micro tube.

SARSTEDT AG & Co. KG

Sarstedtstraße 1
D-51588 Nümbrecht

Phone: +49 2293 305 0

export@sarstedt.com
www.sarstedt.com