









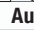
# RF3000<sup>®</sup>

pipette controller

## Packed with features for comfortable and accurate pipetting

The powerful pump offers rapid, yet sensitive aspiration and dispensing, has adjustable valves and is microprocessor controlled for constant reliability and accuracy.

- Operates with 100–240 volt electrical supplies worldwide
- Universal grip and soft coated triggers with direction indicators, for right or left handed users
- Ergonomically balanced in the hand with a loaded pipette
- Bright LEDs indicate three operation modes: high, low, and gravity
- Toggle operation by pressing the large mode button
- Filters are readily available, so you are not tied into a pipette specific brand—helps reduce the total life costs of ownership versus other makes of pipette aids
- An adapter can be fitted to the nose cone to stabilize small volume and thinner pipettes
- Low battery indicator; can be recharged when in use

Part No.	Description
<b>HS3000</b>	RF3000 with power adapter
<b>Replacement parts</b>	
<b>HS5000RFFIL</b>	Filters, pack of 5
<b>HS300501</b>	Universal power adapter
<b>HS300502</b>	Battery RF3000    
<b>HS300504</b>	Nose cones, pack of 5      <b>Autoclavable</b> *Assorted colors (Black, Yellow, Red, Grey and Blue)

Specifications	
<b>Power:</b>	Rechargeable li-Ion battery (replaceable)
<b>Filter:</b>	PTFE/25 mm dia / 0.45 µm filter
<b>Pipettes:</b>	1 – 100 mL, glass and plastic
<b>Pipette Nozzle:</b>	Polypropylene, autoclavable
<b>Dimensions:</b>	15 x 152 x 41 mm
<b>One-year warranty</b>	Tested to comply with U.S., Canadian and European Safety and EMC Standards

